



PRORESET NOVALIGHT

USER MANUAL

PRORESET

Congratulations on choosing one of the most innovative and safe red light therapy beds on the market. Now that you have all the equipment to begin your red light journey, it is extremely important to read this user manual before operating your red light therapy bed.

SAFETY FIRST

At ProReset®, safety is one of our main priorities. It is our policy and practice to ensure that ProReset® red light therapy beds have the highest quality, comfort, and safety.

We have 5 features for your comfort:

- EPISTAR® LED chips with 99.8% light transmittance
- Integrated self-circulating cooling system realizes 360-degree all-round heat dissipation and greatly improves comfort
- Smart control system with built-in timer and status queries
- Automatic session recording for easy review
- Convenient remote access and control from anywhere.

PREFACE

Thank you for purchasing a ProReset® red light therapy bed (hereinafter referred to as "this product")!

This product is constructed from high-quality ABS engineering plastics and aviation-grade aluminum. It also incorporates UK Lucite acrylic sheets and Taiwan EPISTAR® LED chips with 99.8% light transmittance, ensuring durability, safety, and performance. It can be widely used by families in home settings, wellness and community service centers, sports venues, gyms, spas, beauty and health centers, physiotherapy centers, ships, and other public environments.



PRECAUTIONS

To ensure safe use, please read the following precautions before use.

Important Safety Considerations:

Do NOT use the red light therapy bed if you are:

- Pregnant
- Experiencing cold or flu symptoms
- Under the influence of alcohol
- Diagnosed with severe emphysema, malignant tumors, or active tuberculosis
- Suffering from hypertension (above 160/100 mmHg)
- Using a pacemaker without prior approval from the manufacturer
- Taking photosensitizing medications or experiencing photosensitivity
- Experiencing active hemorrhage, organ damage, or severe cardiovascular/kidney disease.

Warning:

- If an abnormal smell occurs during use, turn off the power and pull out the plug
- Please make sure that the socket plug is fully inserted firmly
- Please do not disassemble the product at will
- This product should be kept away from direct exposure to sun, water and fire sources
- Children must be under the supervision of adults at all times.
- Please do not place heavy objects on top of the red-light therapy bed or any accessory
- Do not use if any wire is loose
- Do not pull or twist the power cord forcibly
- Do not touch or plug the power plug into an electrical outlet with wet hands
- Do not use the red-light therapy bed inside a bathroom, near a pool or other humid environments.

PRECAUTIONS FOR USE

1. It is recommended to irradiate naked, and all accessories and cosmetics must be removed from the whole body.
2. Pay attention to eye protection during use. Please remove contact lenses during use, and special eye protection (eye patch) must be worn.
3. The use time should be set by a special person according to the suggestion of the physical therapist. Do not adjust the irradiation time at will during use.
4. Pay attention to the user's reaction during use, such as overheating, palpitation, dizziness, etc., and contact the staff immediately.
5. If there is a wound, the wound should be debrided before use, and also used with caution on fresh scars.
6. After the first use, observe whether there is any red light allergic reaction. If not, continue to follow the course of treatment.
7. It cannot be used to treat congenital skin hypoplasia, such as warts, spots, and moles.
8. Under normal circumstances, each use should not exceed 30 minutes.

SPECIFICATION STRUCTURE

1. Control panel
2. Red light therapy bed body
3. Chassis



PRODUCT INSTALLATION

Unpacking Inspection

Please check if the packaging is intact before unpacking. After unpacking, take out the parts, please check the following items:

- 1 Confirm that there is no damage to the product during transportation, and check that the appearance of the product is not damaged or notched;
- 2 Check whether the product nameplate information is consistent with the ordered product;
- 3 If you have ordered accessories, confirm that the received accessories are what you need.



WARNING:

Do not install or operate any photodynamic therapy bed that is damaged or has faulty parts, otherwise there is a danger of injury. If you find any questions, please get in touch with ProReset® and solve them immediately.

Product Installation Specifications

- 1 The product must be grounded correctly before using. When multiple products are grounded, it is recommended that you do not use a common ground wire to avoid the ground wire forming a loop;
- 2 In order to prolong the service life of the lamp and the product, please turn off the power in time when not in use;
- 3 After the main power switch is turned off, you must wait at least 3 minutes for the product to cool down completely before allowing installation/ disassembly operations;
- 4 Never install or test the unit without proper insulation, as this may result in electric shock;
- 5 The indoor air must be kept unblocked, and the distance between the outer edge of the product and the wall is not less than 30cm;

- 6

Please review the installation, usage and environment requirements outlined in this manual. Failure to comply may shorten the product's service life or result in damage;
- 7

This product is delivered as a completed device. After opening the package, it can be used after wiring according to this manual. Do Not disassemble or repair the product without electrician or seller's guidance.

Specification structure

Material	ABS engineering plastics and aviation aluminum, UK Lucite acrylic sheet
Size (LxWxH)	227.5 x 123.5 x 112 cm
Weight	470 kg
Light Source	China Taiwan EPISTAR® LED chips
Led Quantity	41,600 pcs LEDs
Irradiance	188 mw/cm²
Total Power	5,200 W
LED Lifetime	100,000 hours
Wavelength	633 nm, 660 nm, 810 nm, 850 nm, 940 nm
Voltage	220 V / 380 V
Frequency	1-10,000 Hz
Control System	Smart Control System 2.0, (Optional External Controller)
Features	Self-circulating air-cooling system, surround sound Bluetooth audio system

INSTALLATION REQUIREMENTS

This product needs to be installed indoors. During installation, there is no obstruction within 1 meter around the product to ensure well heat dissipation.

Temperature and Humidity

- The average temperature must be between 15°C and 35°C;
- The relative humidity of the air is less than or equal to 75%RH, without condensation.

Other Installation Environment Requirements

- Install in a place that is unlikely to be subjected to severe vibration and impact;
- Install in a place where metal powder, dust, oil, water, etc. cannot penetrate into the product;
- Avoid installing in places with oil mist, steam, or salt;
- Avoid installing in places directly exposed to the sun;
- Avoid installing in dusty or metal powder places;
- It is strictly forbidden to install in places with flammable and explosive materials such as anesthetics.



NOTE:

If there are special installation environment requirements, please consult and confirm in advance.

Installation Space and Distance Requirements

- The min distance between the four sides of a product and the wall – 30cm.
- The minimal distance between two products – 100cm.

Electrical Configuration Requirements

This product MUST BE equipped with a leakage protection switch or an air switch; otherwise there is a danger of injury. This device is not provided by our company, please configure by yourself after consulting. The power cord of the device must be directly connected to the junction box of the dedicated circuit (it is HIGHLY RECOMMENDED to be connected by a professional electrician, refer to the power configuration of the local and/or country code, and be equipped with appropriate circuit breakers and wire gauges).

If the power supply is configured with 220V:

It is necessary to configure a dedicated single-phase circuit capable of providing 40 amp at the installation location.

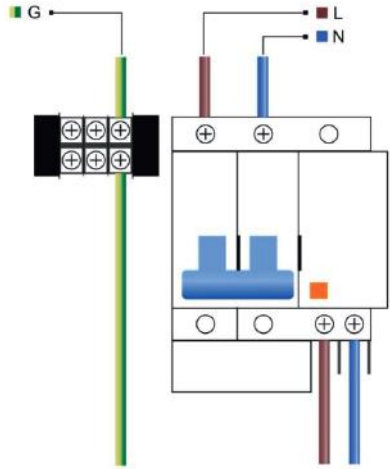


Figure 1: 220V, 40 amp

If the power supply is configured with 380V:

It is necessary to configure a dedicated three-phase circuit capable of providing 3*25 amp at the installation location.

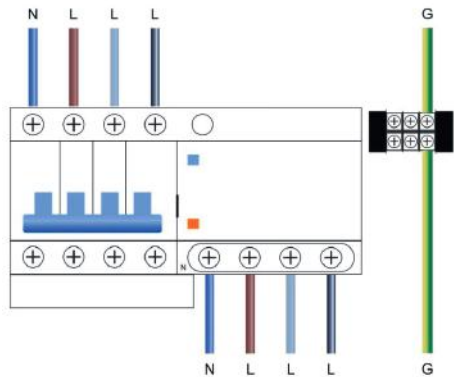


Figure 2: 380V, 3*25 amp

PRODUCT USE AND OPERATION

Please read the operating instructions carefully before using the product, and operate in strict accordance with the operating instructions.

Internal control system operation panel

Main Screen

After the device is powered on, enter the main screen, as shown in the picture.

Each numbered component on the screen serves the following function:



1-4 buttons

Working time adjustment button (It cannot be adjusted after starting work.)

- 1 button** press to increase time by 1 minute.
- 2 button** press to decrease time by 1 minute.
- 3 button** press increase time by 1 second.
- 4 button** press to decrease time by 1 second.

5 – Session time display

It will display the currently set working time when adjusting. When working, it will display the remaining working time (countdown).

6 – The current fan gear display

When the fan is off, the icon does not move, and there is no digital display in the middle. When the fan is turned on, the icon rotates, and the current gear of the fan is displayed in the middle of the icon.

7 – Equipment operation control button

Press "START" icon to start the session, press the "PAUSE" icon to pause the session, and press the "CONTINUE" icon to continue the session. After starting the session, long press the button for 2 seconds to end the work. It can be operated freely when it is not connected to the wireless controller. After connecting to the wireless controller, it is necessary to click "START" icon on the home screen of the wireless controller to allow the operation.

8 – Setting and Music Control button


Press to enter the music control screen. Long press for 3 seconds to enter the setting screen.

The icon is displayed in grey in the middle when disconnecting the Bluetooth and turns to blue when connected.

Music play control screen

The functions of each digital marking part on the screen are as follows:

1. Play / Pause button
2. Previous song
3. Next song
4. Volume Down
5. Volume Up
6. Current Volume
7. Back to main screen
8. WiFi Connect Icon

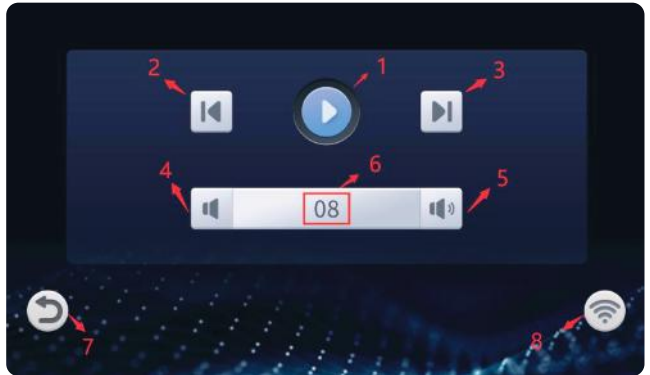

The wireless connection icon has the following seven states:



Not connected to the network.



Indicates that the network is connected but the server is not connected and the signal is poor.



Indicates that the network is connected but the server is not connected and the signal doesn't fully reach the product.



Indicates that the network is connected but the server is not connected and the signal is good.



Indicates that the network and the server are connected. The signal is poor.



Indicates that the network and the server are connected. The signal doesn't fully reach the product.



Indicates that the network and the server are connected. The signal is good.

Pulse setting screen

As shown in the picture, this screen can configure the pulse output of the device.



You can adjust the Frequency and Duty settings by putting values in the corresponding boxes and pick output status by checking the box for each channel individually. Alternatively, you can select multiple outputs, enter the desired values in the red box above, and press the SET button to apply the specified Frequency and Duty settings to all selected outputs simultaneously.

WiFi Setting Screen

This screen has two functions:

a) WiFi configuration:

Click the "Net name" edit box, and on the screen that pops up, select the name of the network to be connected, and enter the corresponding password in the "Password" edit box, and press "START" to start reconnecting the WiFi network.



The system will save the corresponding passwords of the 9 most recently connected network names. If the password edit box is not empty after selecting the network, it means that the network has already been connected. At this time, there is no need to enter the password again, and you can directly press the START button to connect the WiFi network.

b) Local configuration:

You can modify the device name in the "Device name" edit box, and modify the local UDP broadcast receiving port in the Local port edit box. The port must be consistent with the UDP broadcast port set on the APP to connect to the corresponding APP.

Usage records Screen

As shown in the picture above, you can query the latest 999 usage records on this screen.

Press "Last" to query the previous time, and press "Next" to query the next time.



Date setting screen

After filling in the correct date and time on this screen, press the "OK" button to set the information, or directly press the return button to return to the main screen without saving the settings.



Fan setting screen

Each numbered component on the screen serves the following function:



1. Set the gear of the cooling fan after the session.
2. Set the temperature at which the product's cooling system starts to work.
3. Set the product's temperature at which the product overheat warning system appears.
4. Set the delay time at which the cooling fan turns off after a session.

WARNINGS.

IMPORTANT READ BEFORE USE.

Safety Requirements

- The safety of the product can only be guaranteed when it is connected with the product provided by the company, and the safety of the equipment provided by any other company cannot be guaranteed;
- In order to prevent circuit danger, the product can only be connected to a power socket with protective grounding. If the power socket is not connected to a grounding wire, or the integrity of the protective grounding wire is questioned, please do not use the socket;
- When a sudden power failure occurs during the normal use of the product, the product will not cause damage to the operator or user. Please turn off the power switch first, and perform normal operation according to the instructions when the power is turned on again;
- Do not open the casing of the product by yourself; otherwise, there may be a risk of electric shock. The maintenance or upgrade of the product must be carried out by maintenance personnel trained or authorized by the company;
- To prevent electric shock and reduce the risk of product failure, do not expose the product to water.
- If the product is accidentally exposed to water, please stop using it immediately and hand it over to professional technicians before using it again;
- Before using the product, check the product and its connecting cables and accessories carefully to ensure that they can work normally and safely;
- Use according to the parameters specified by our company;
- In order to ensure the safe operation of the product, all kinds of replaceable parts, accessories and various consumables equipped with the product please use the products provided or designated by our company;

**If you have any questions, please contact us at
support@proreset.com**

Disclaimer

ProReset® NovaLight is not intended to diagnose, cure, or prevent specific diseases or medical conditions. These products are only intended to be used for relaxation and wellness purposes.

Any information contained in this manual is not intended as medical advice or a substitute for medical advice.

ProReset® does not intend these products to cure any physical or skin conditions solely from their use. A qualified health professional or doctor should always be consulted prior to using the products, and professional medical attention should be sought regarding any serious or life-threatening conditions.

These statements have not been evaluated by any Food and Drug Administration or health authorities. This product is not intended to diagnose, treat, cure, or prevent any disease whatsoever.

ProReset® is not responsible for any damage or loss that may arise from the use of our products, services, or information, and we disclaim all liability to the fullest extent permitted by law.

Maintenance and Storage

- Clean the cabin surface with alcohol once a week to maintain hygiene and performance
- Store the device only when it has cooled down to room temperature.
- Keep it in a dry environment, away from direct sunlight.
- Protective covers should be used if available.
- Do not store near flammable items.
- Wipe the surface with a damp cloth.
- Follow local disposal practices as they may damage your device.
- Check that the fan is running normally and no debris is stuck, etc.
- Check whether the input power voltage complies with the standards.
- When not in use, keep the door open and unplug the power cord
- When maintaining or cleaning the interior of the operation cabin, ensure that the power supply is disconnected, and do not use liquid to clean the product parts;

- If the wire is damaged, it must be replaced by the manufacturer/service agent/related professionals to avoid injury;
- The organic glass plate in the engine room is damaged, cracked or scratched, please replace it immediately;
- The storage must be placed in a dust-free, dry environment;
- The ambient temperature of the storage location must be within the range of 15°C~35°C;
- The relative humidity of the storage location must be 0%~75%RH, and there should be no condensation;
- The storage location should avoid corrosive gases and liquids;
- There should be proper packaging protection, and should be stored above the ground.

Disposal

This device must not be disposed of with your normal household rubbish. This device requires disposal via your designated waste collection point.



ProReset® Warranty Policy

ProReset® warrants that its products are free from defects in materials and workmanship under normal and intended use. This 3-year limited warranty begins from the date of production and covers the original purchaser only. This warranty is non-transferable and requires proof of purchase for any claims. Annual maintenance, inspection and/or calibration services must be performed by ProReset® or an authorized service provider to maintain warranty validity. Failure to comply will void this warranty.

Warranty Coverage

ProReset® will, at its sole discretion, repair or replace defective products or components found to be faulty within the warranty period.

Parts not manufactured by ProReset® but included in the product carry the original manufacturer's warranty.

EXCLUSIONS

This limited warranty does not cover:

1. Normal wear and tear.
2. Consumable parts, including air filters.
3. Shipping costs to and from ProReset's warehouse for repairs or replacements.

4. Damage incurred during transit (claims should be filed with the carrier).
5. Accessories or components added after shipment that are not approved by ProReset®.
6. Products purchased through unauthorized dealers or third-party resellers.
7. External accessories not listed in ProReset's standard product catalog.
8. Yearly or other maintenance or service costs beyond warranty repairs.

WARRANTY VOIDING CONDITIONS

This warranty will be void if:

1. Cosmetic damage occurs not related to material or workmanship defects (e.g., scratches, dents, or tears).
2. Damage results from accidents, misuse, abuse, improper use or storage, or unauthorized modifications.
3. Repairs or alterations are performed without ProReset's prior written approval.
4. The product is used in conditions exceeding the guidelines outlined in the ProReset Operating Manual.
5. The serial number is removed, defaced or untraceable.
6. The product is used with non-approved accessories.
7. The product contains any chemicals (e.g., exposed paint thinner) during use.
8. The product is placed too close to heat sources, exposed to direct sunlight or solvents.
9. The product is not operated within the specified operational conditions.

Claim Process

- Claims must be submitted within 30 days of defect discovery and no later than the warranty expiration.
- Claims must include a detailed description of the defect with proof of purchase.
- Purchasers are responsible for shipping costs, import costs, and return to/from ProReset's warehouse for warranty repairs.
- Claims must be submitted via certified mail or through the ProReset website.

Limitations of Liability

- ProReset's liability is limited to repair or replacement at its discretion under this warranty.
- No compensation is provided for incidental damages, loss of use, downtime, or any consequential losses.

General Conditions

- Purchasers assume full responsibility for compliance with local safety codes and operational regulations.
- This warranty applies only to products used and maintained per ProReset's instructions and specifications.
- Resale, rental, lease, or transfer of ownership voids this warranty.
- No implied warranties of merchantability or fitness for a particular purpose extend beyond this document.

Troubleshooting

The following issues may be encountered during use, please refer to the following methods for simple troubleshooting analysis:

Display doesn't turn on

Use a multimeter to check whether the input power voltage of the product is normal.

The power air switch trips after power-on

Check whether there is a grounding or short circuit between the input power supplies, and eliminate the problem.

The lamp on the display does not work

- Check whether the interface tube is loose;
- Check whether the tube is burned out, if so, please seek service.

The content in this manual is subject to change.

Download the latest version by scanning this QR code
or from <https://proreset.com/manuals-xj37tyyq-proreset/>



**If you have any questions about your product or this user manual,
please contact ProReset® by sending an email to
support@proreset.com.**



ProReset® S.L.
C/Londres 18 - P.I. Mas de Taus
La Pobla de Vallbona Valencia 46185 Spain
WhatsApp: +34 624 198 129
Phone: +34 865 699 306
support@proreset.com
www.proreset.com